

DSC & Modulated DSC (MDSC®) Theory and Applications Online Courses

Part 1: Theory & Instrumentation



©2020 TA Instruments

1

Differential Scanning Calorimetry (DSC)



2

Agenda

- Understanding DSC and MDSC techniques
 - How heat flow is measured on a DSC
 - When and why it is important to use MDSC
 - How MDSC heat flow signals are calculated
- Calibration & Verification
- Instrument & Method Considerations
 - Purge Gas
 - Cooling Accessories
 - Sample Pans & Sample Preparation
 - DSC Cell Cleaning
 - Experimental Setup and Method Development


3



3

What is a Differential Scanning Calorimetry

- A DSC measures the difference in Heat Flow Rate between a sample and inert reference as a function of time and temperature

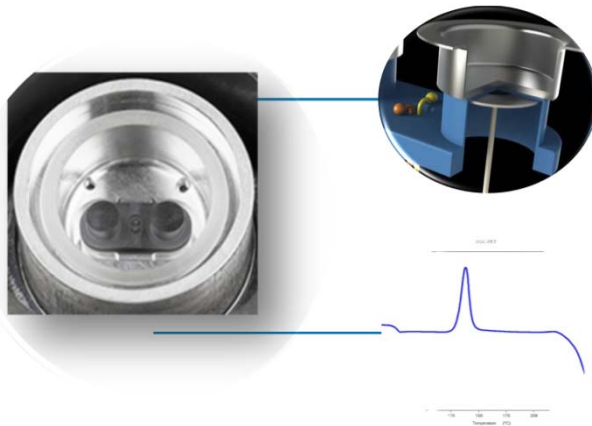


4



4

Simple Heat Flux DSC Cell Schematic

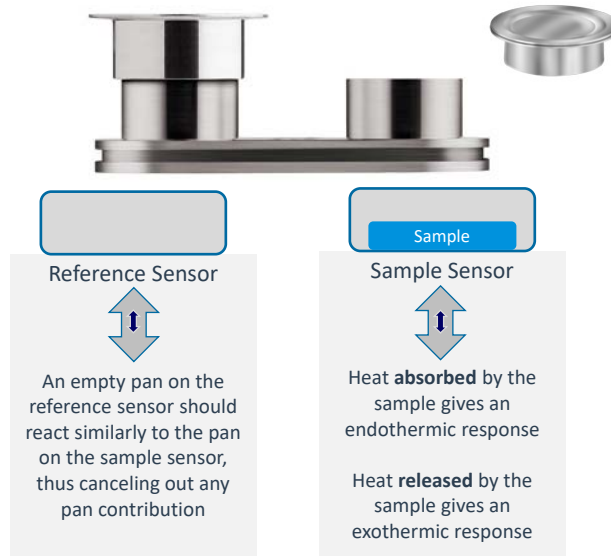


5



5

Endothermic Heat Flow – Heat Absorbed by Sample

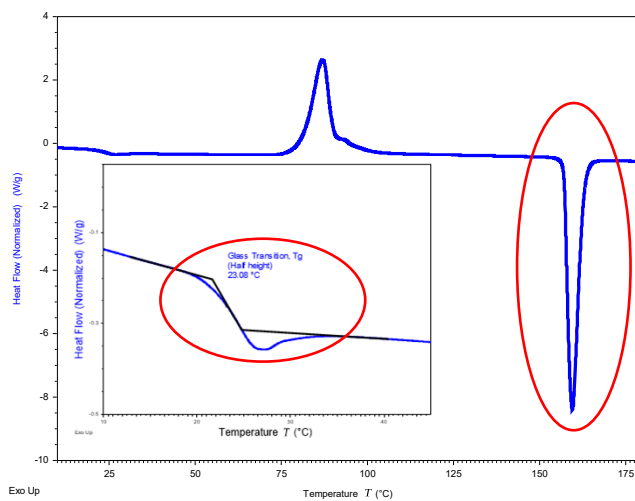

Endothermic Events


- Glass transition
- Melting
- Evaporation/volatilization
- Enthalpic recovery
- Polymorphic transitions
- Some decompositions

6



6

Exothermic Heat Flow – Heat Released by Sample

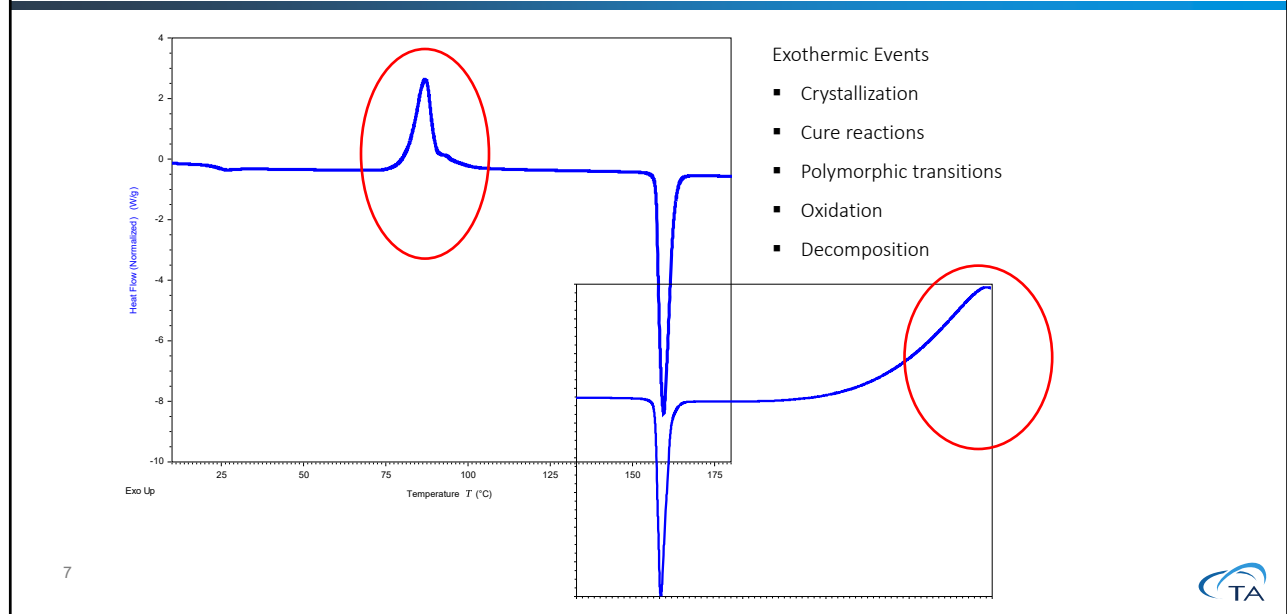

7

DSC Heat Flow Equation

$$\frac{dH}{dt} = \text{DSC heat flow signal}$$

(mW or $\frac{\text{mJ}}{\text{s}}$)

$$C_p = \text{Sample Heat Capacity}$$

= Sample Specific Heat x Sample Weight

$$\frac{dH}{dt} = C_p \frac{dT}{dt} + f(T, t)$$

$$\frac{dT}{dt} = \text{Heating Rate } (^\circ\text{C}/\text{min})$$

$$f(T, t) = \text{Heat flow that is function of time}$$

at an absolute temperature (kinetic)

8



8

DSC Heat Flow

$$\frac{dH}{dt} = C_p \frac{dT}{dt} + f(T, t)$$

Heat
Capacity

Glass Transition
Some Melting

Kinetic

Crystallization
Some melting
Cure reactions
Volatilization
Decomposition
Denaturation

9



9

What DSC Can Tell You

- Identification of amorphous & crystalline material
- Identification of phase transitions and changes of state
- Specific Heat Capacity – measure of molecular mobility
- Heats of fusion & reactions
- Oxidative/Thermal Stability
- Reaction Kinetics

10



10

Modulated DSC[®] Theory (MDSC[®])



11

What is MDSC?

- MDSC separates the Total heat flow of DSC into two parts based on the response of the system to a changing heating rate. The changing heating rate is effected by 'superimposing' a sinusoidal heating rate on a linear heating rate.
- In general, only heat capacity and melting respond to the changing heating rate resulting in an increase in signal. Kinetic events tend to occur at different temperatures as a function of heating rate. For example, increasing the heating rate can shift decomposition to a higher temperature.
- The Reversing and Nonreversing signals of MDSC are not necessarily a measure of reversible and nonreversible properties.

12



12

MDSC[®] Theory: Heat Flow Signals

$$\frac{dQ}{dt} = C_p \frac{dT}{dt} + f(T, t)$$

Total
Heat Flow

- All Transitions

Reversing Heat
Flow

- Heat Capacity
- Glass Transition
- Melting

Non-Reversing
Heat Flow


- Enthalpic Recovery
- Evaporation
- Crystallization
- Thermoset Cure
- Denaturation
- Decomposition
- Some Melting
- Chemical Reactions

13



13

Isothermal @ 25°C

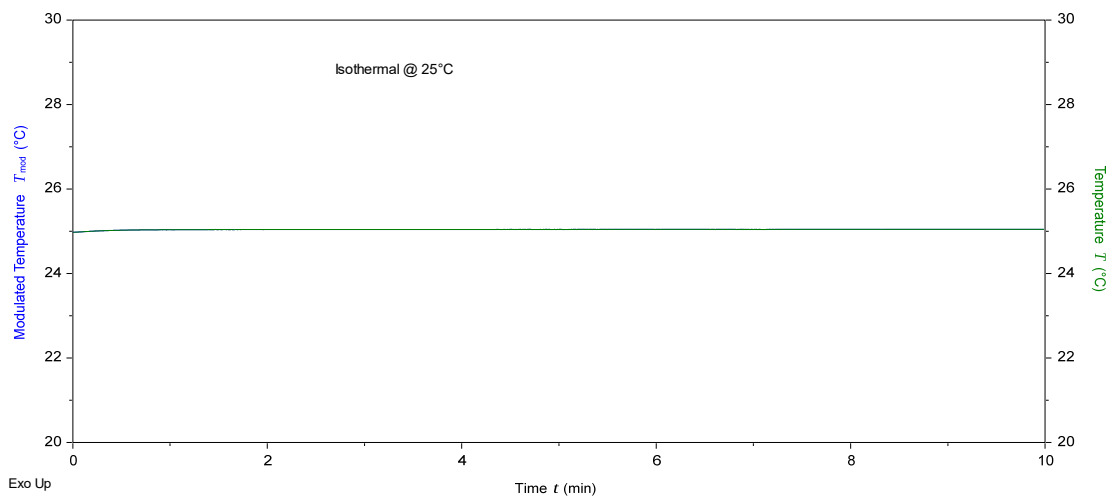


14




14

MDSC - Quasi-Isothermal @ 25°C

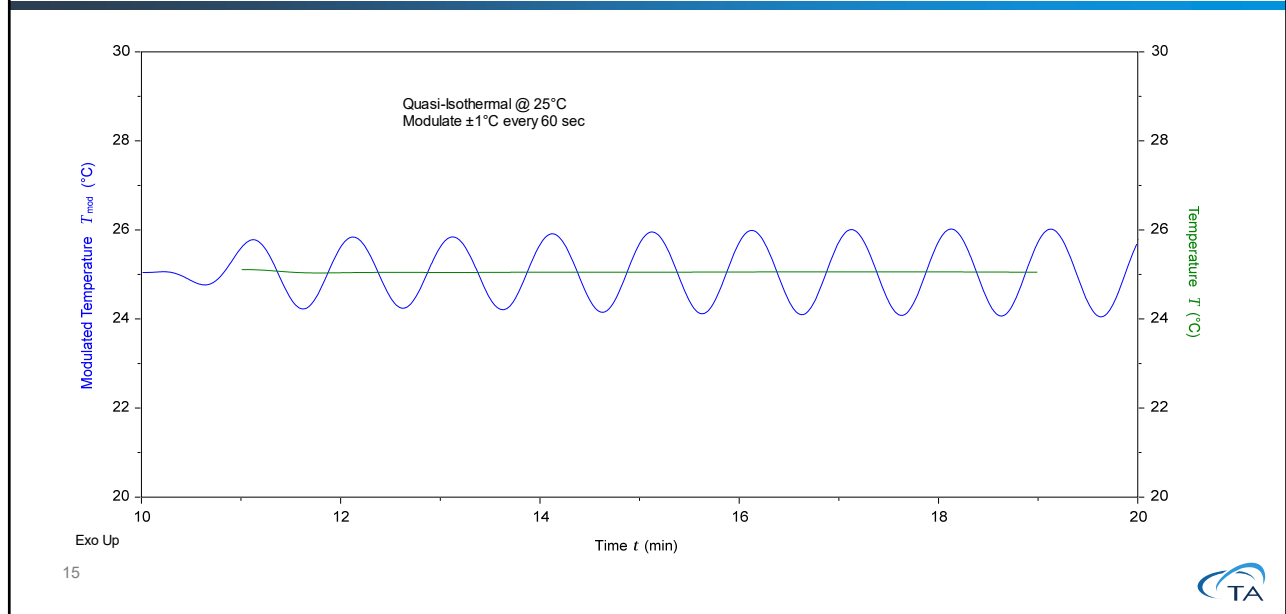

15

MDSC – 1°C amplitude, every 60 sec, @ 2°C/min

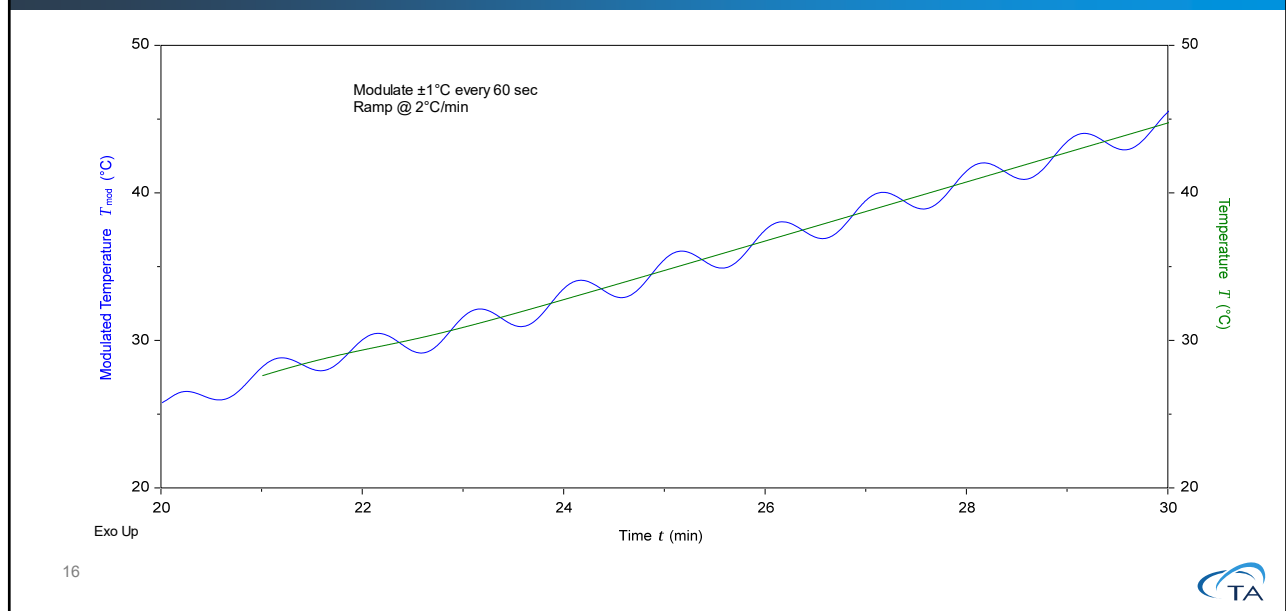

16

MDSC[®] Theory: Calculation of MDSC[®] Signals


- All MDSC signals are calculated from three measured signals.
 - Time
 - Modulated Temperature and by implication Modulated Heating Rate
 - Modulated Heat Flow
- During the setup of the MDSC experiment, the user enters the following parameters:
 - Average (or underlying) heating rate (°C/min)
 - Temperature modulation period (seconds)
 - Temperature modulation amplitude (°C)

17



17

Average & Modulated Temperature

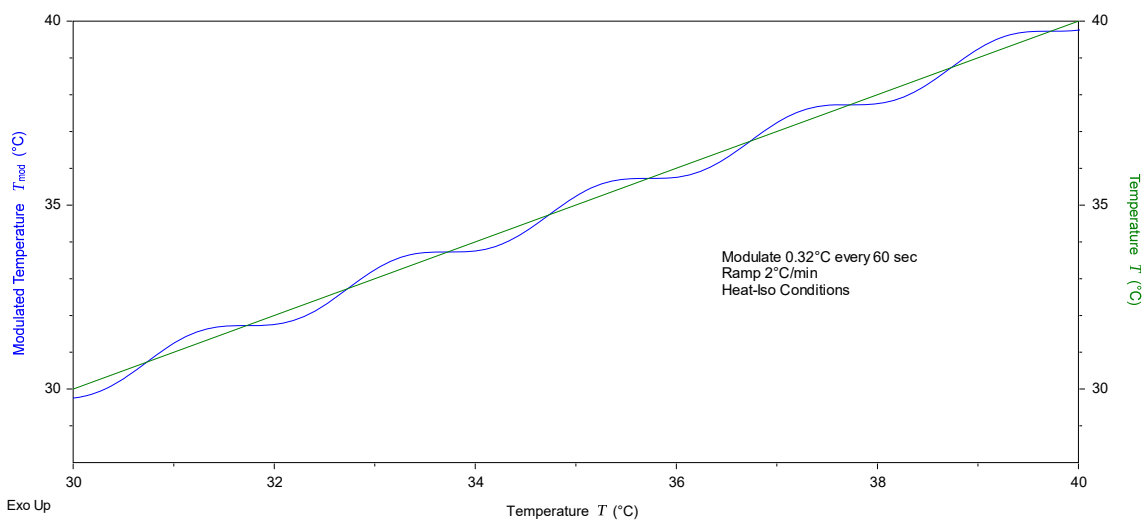


18



18

Average and Modulated Heating Rate

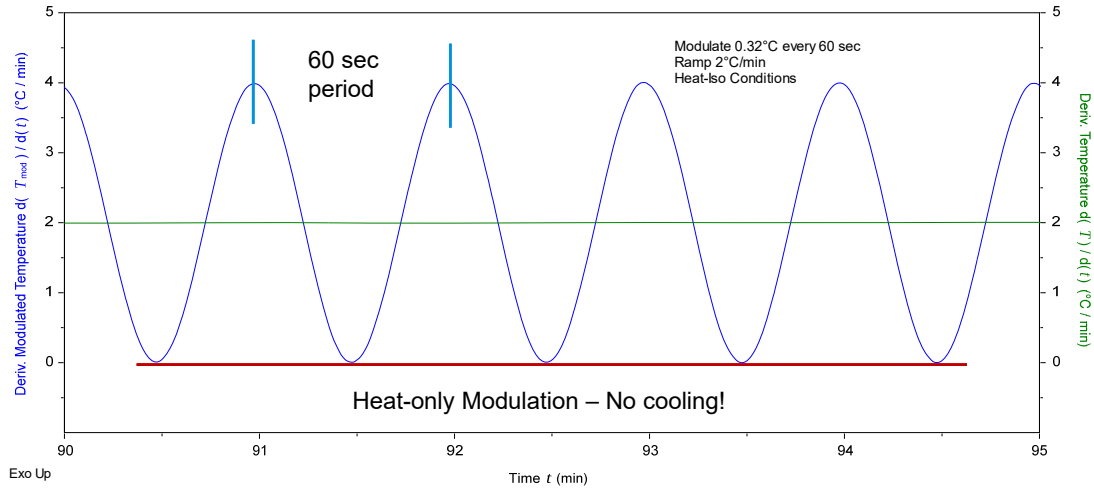


19



19

MDSC® Raw Data Signals

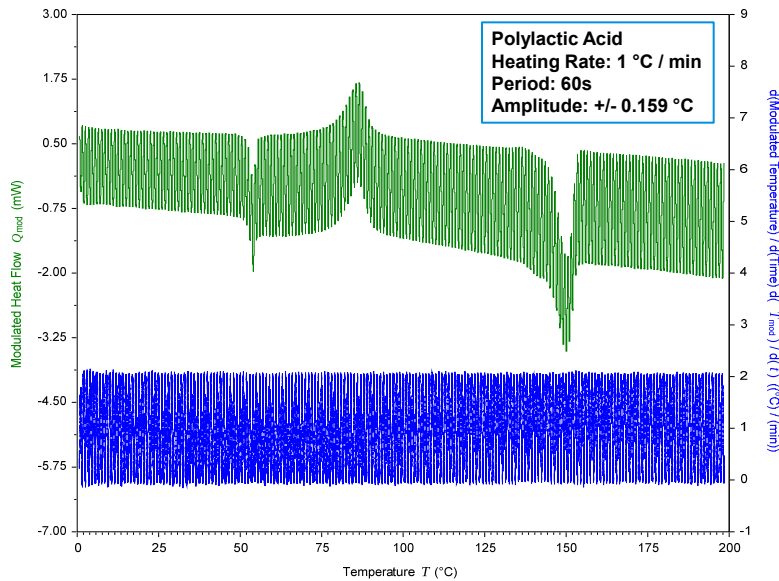

20



20

MDSC Theory: Calculation of MDSC® Signals – Simple Deconvolution


- Raw data is averaged over a period of 1 oscillation and the average is subtracted from the raw data.
- Modulation is analyzed using a Fourier Transform which yields the amplitude of the heat flow response at the modulation frequency.
- This results in the following:
 - $\langle dQ/dt \rangle$ = average heat flow; Q = heat
 - A_{HF} = Amplitude of heat flow modulation
 - A_{HR} = Amplitude of modulated heating rate

21



21

MDSC® Raw Data Signals: Modulated Heat Flow and Modulated Heating Rate with Calculated Total (or Average) Heat Flow and Average Heating Rate

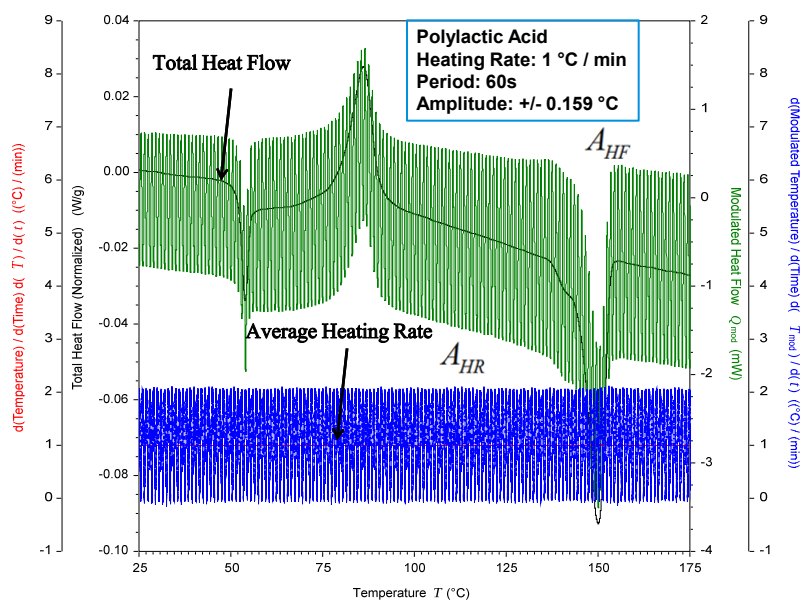

22



22

Reversing Cp Signal

- The reversing heat capacity can be calculated as:

$$Cp_{Rev} = \frac{A_{HF}}{A_{HR}}$$

- Reversing Heat Flow = $Cp_{Rev}\beta$

– Where

Cp_{Rev} = the reversing heat capacity

A_{HF} = Amplitude of the modulated heat flow

A_{HR} = Amplitude of the modulated heating rate

β = Average heating rate

23



23

MDSC Theory: Heat Flow Signals

$\langle dQ / dt \rangle = \text{Average Heat Flow}$

$Cp_{Rev}\beta = \text{Reversing Heat Flow}$


$\langle dQ / dt \rangle - Cp_{Rev}\beta = \text{Non - Reversing Heat Flow}$

24




24

MDSC® Theory: Calculation of MDSC® Reversing Heat Flow and Reversing Cp

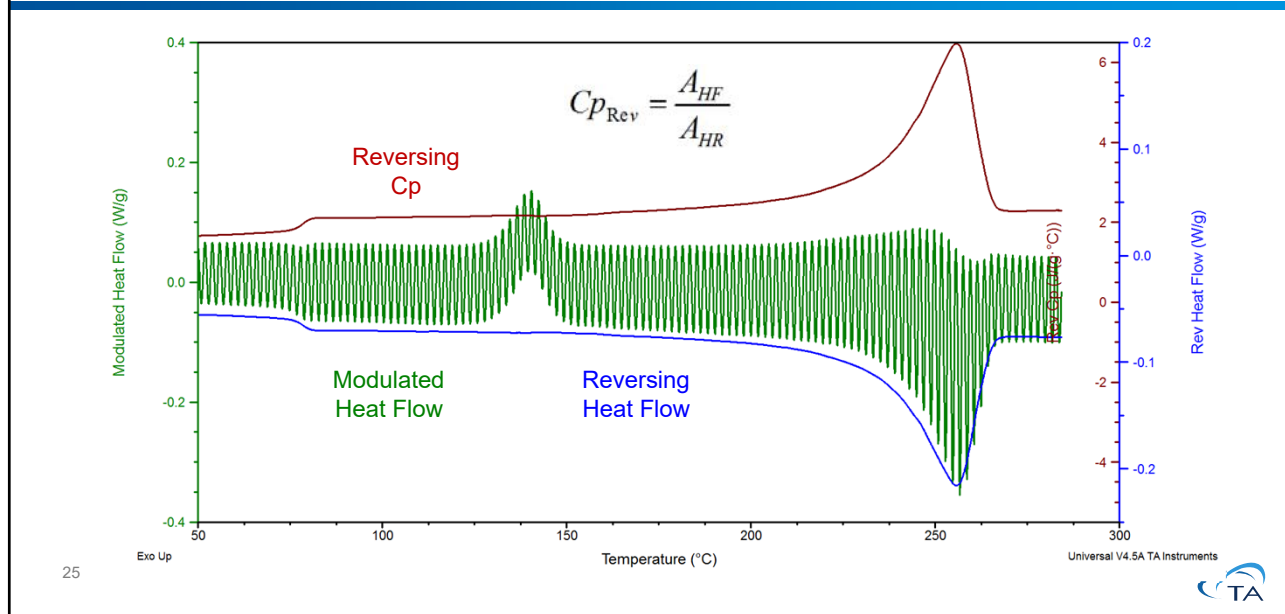

25

MDSC® Theory: Calculated MDSC® Heat Flow Signals Summary

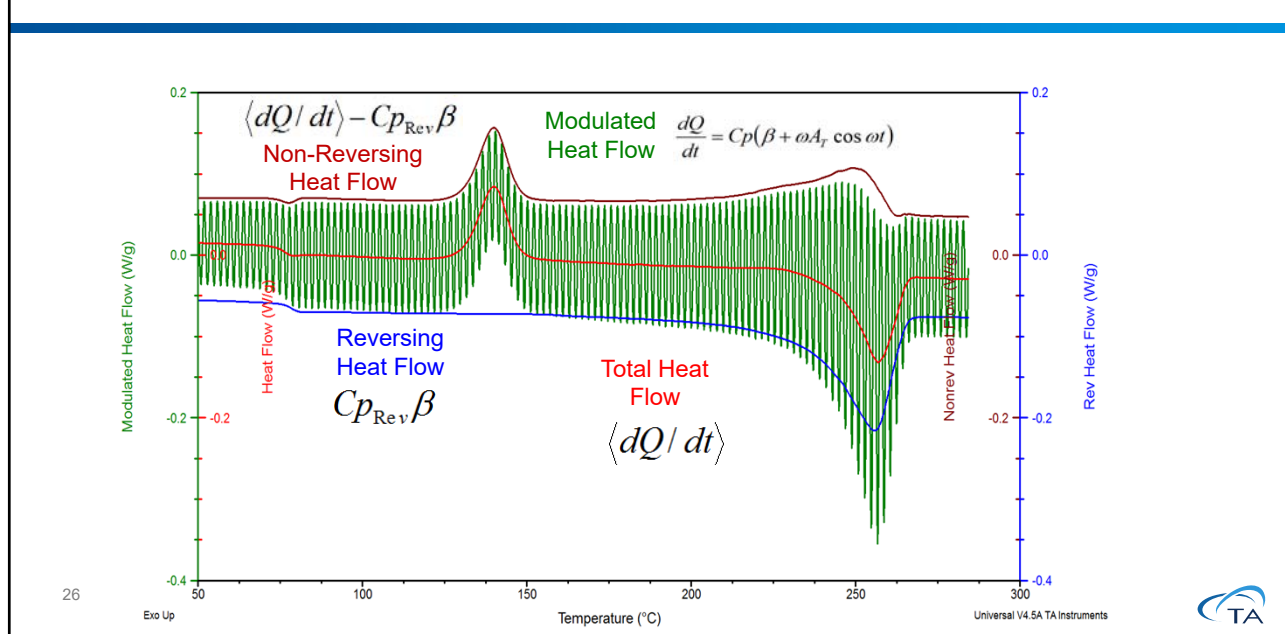

26

Modulated DSC[®] (MDSC[®])

When and Why to Use MDSC



27

When & Why to Use MDSC?

- Run a conventional DSC experiment @ 10°C/min first
 - It may provide all the information you need
- Reasons to run MDSC:
 - Identify heat capacity and kinetic transitions
 - Separate overlapping thermal transitions
 - Detect weak glass transitions
 - Most accurate determination of polymer crystallinity
 - Accurate measurement of heat capacity in a single experiment
 - Gain insight into structural change
 - And many others...

28



28

When & Why to Run MDSC[®] – Amorphous Materials


- If you're looking for a T_g –
 - If the T_g is detectable and can be routinely analyzed, then you may not need to use MDSC
 - However, if the T_g is hard to detect, or has a large enthalpic recovery, then run MDSC

29



29

Epoxy Cured 48 Hours: Heat Cool Heat

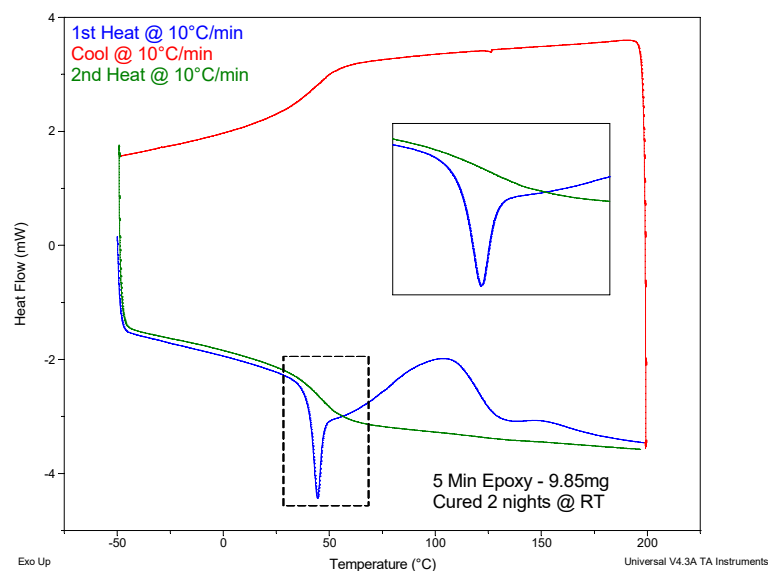


30



30

Rev-Heat Flow Easily Shows Tg

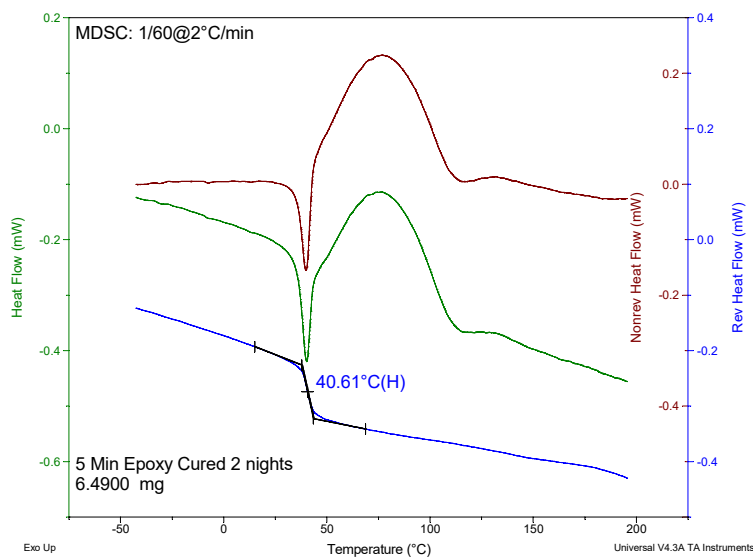

31



31

When & Why to Run MDSC® : Crystalline Materials


- If you are studying polymer melting and crystallization
 - If the melting process appears normal (single endothermic peak) and there is no *apparent* crystallization of the sample as it is heated, then there *may* be no need to use MDSC
 - However, if melt is not straightforward, or crystallization may be occurring as the sample is heated, use MDSC

32




32

Calculated MDSC Heat Flow Signals

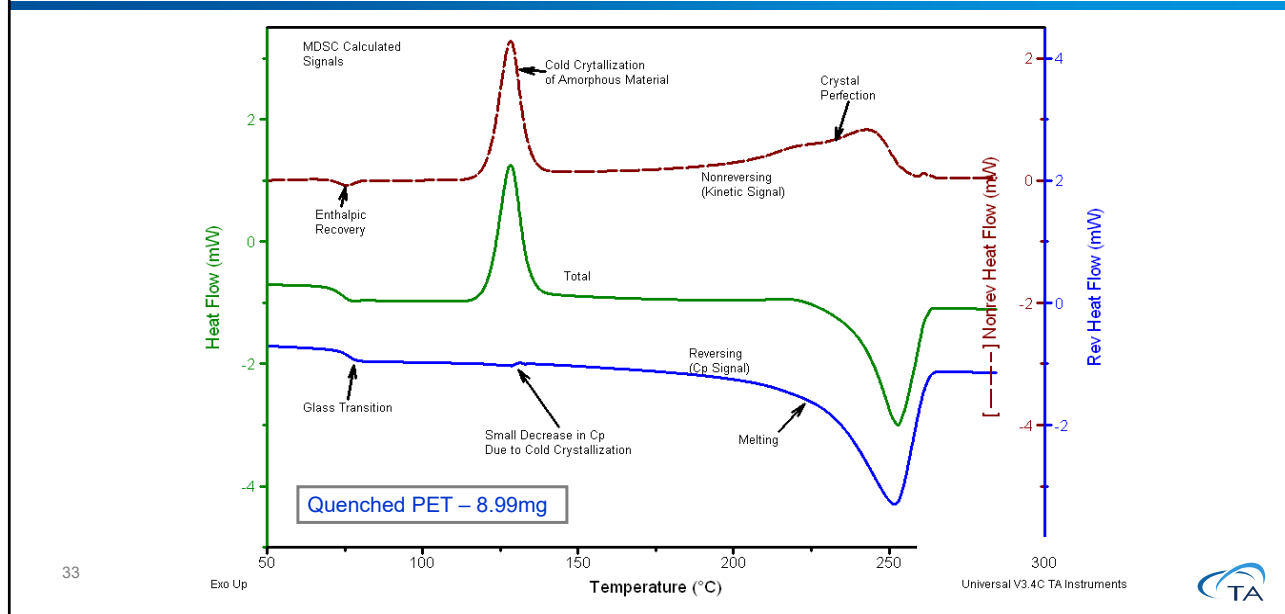


33

DSC of Complex Polymer Blend

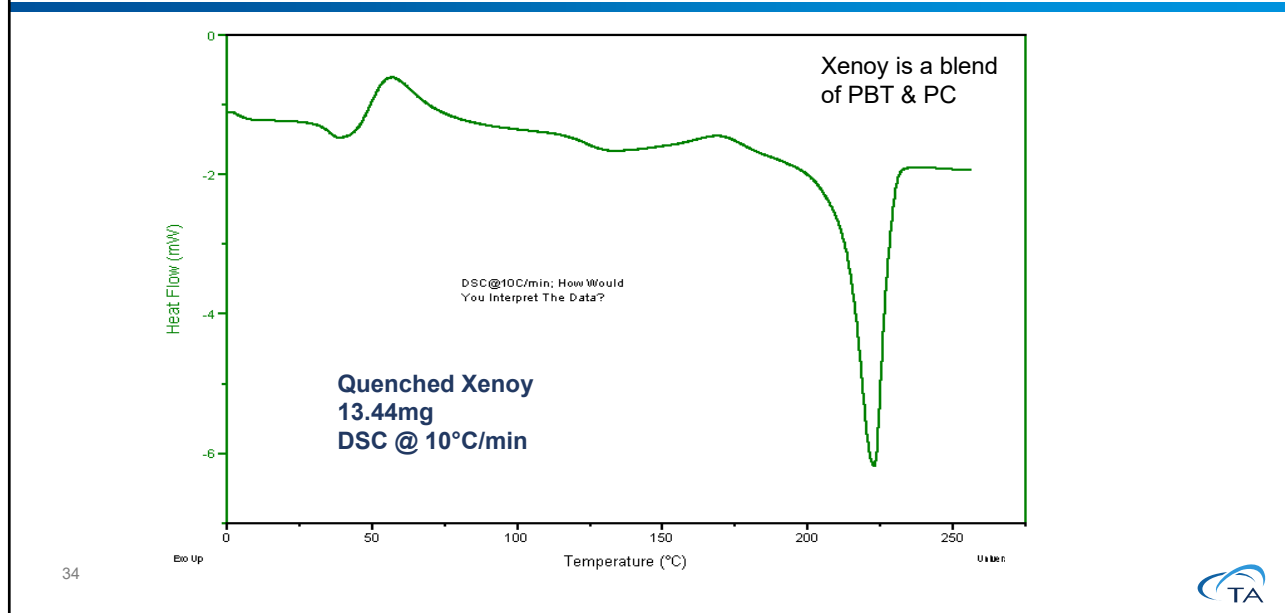

34

MDSC[®] of Complex Polymer Blend

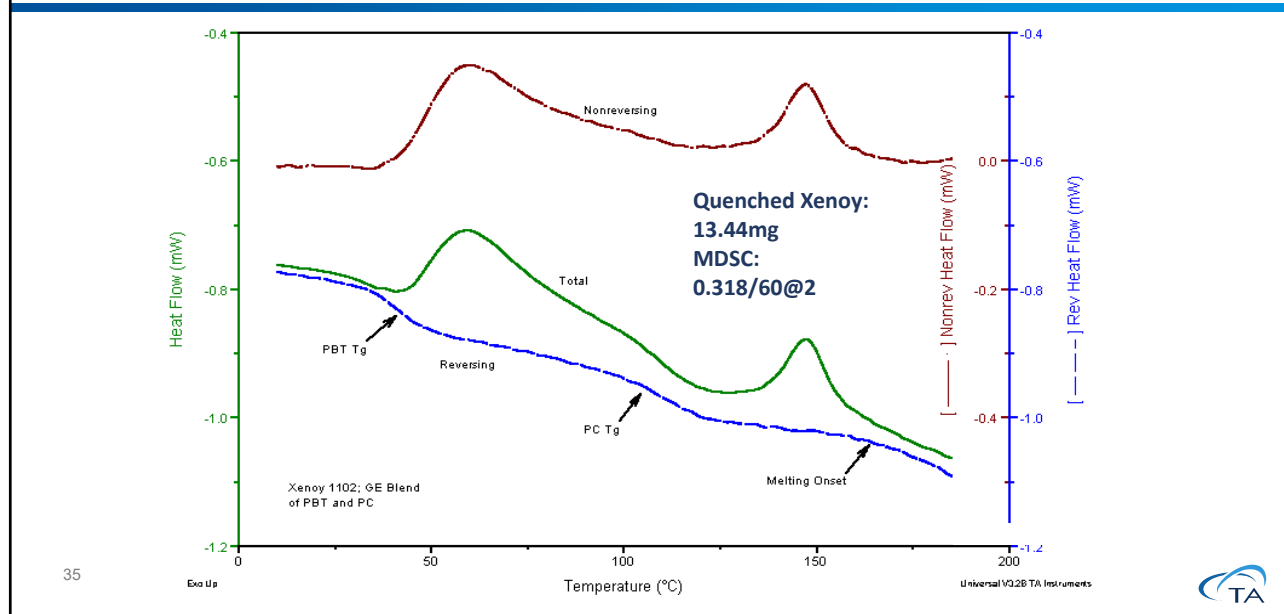

35

When & Why to Run MDSC[®]: Heat Capacity measurements


- If you want to measure accurate heat capacity (C_p), or the change in C_p as a function of time at an isothermal temperature – run MDSC
- Modulated DSC offers two methods for obtaining heat capacity
 - Dynamic Ramp: Use the 'Reversing Heat Capacity' signal, after obtaining a K_{Cp} value.
 - Quasi Isothermal: Very simple and accurate means for measuring specific heat capacity. No calibration is needed, K_{Cp} is determined at each temperature so accuracy is excellent.

36



36

Most Accurate Cp Determination

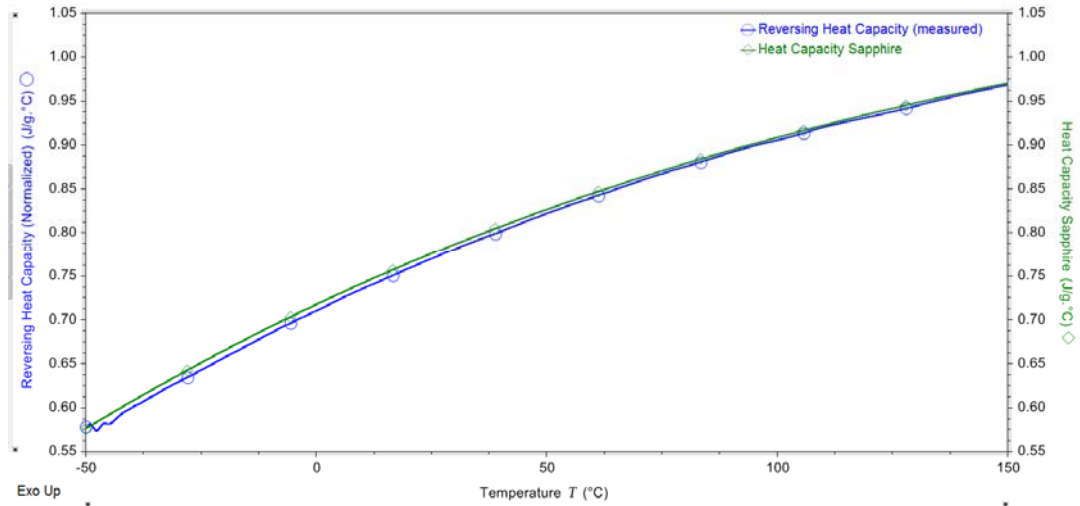

37




37

Calibration & Verification



38

Calibration of Specific Instrument Models

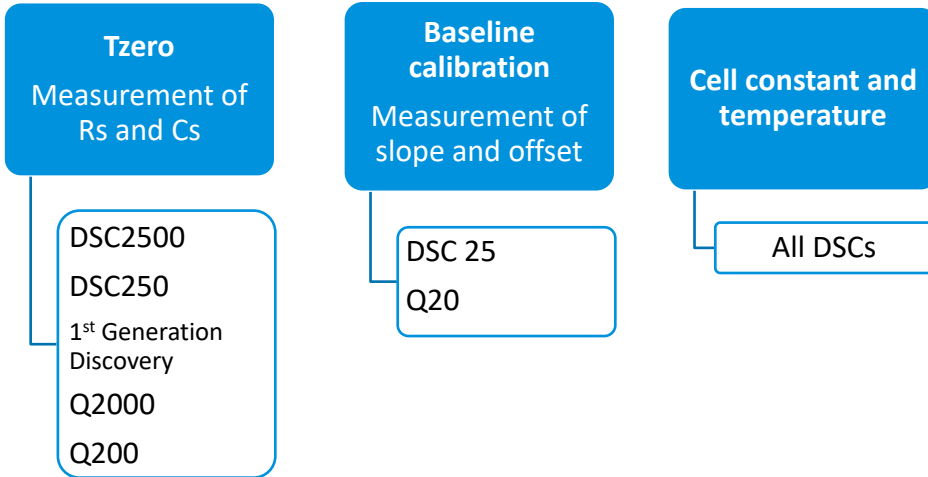


39



39

Calibration Setup in Trios

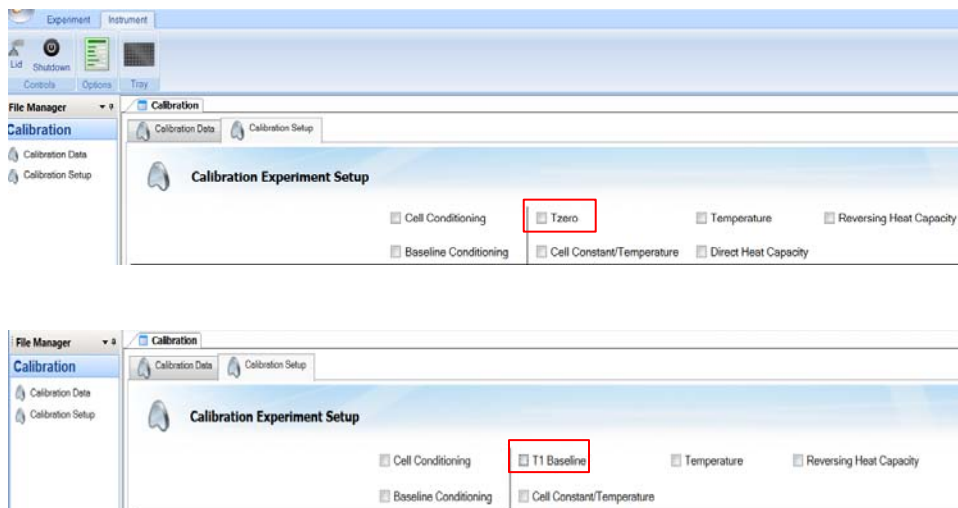

40



40

General Calibration and Verification Guidelines

Calibration

- Use Calibration Mode
- Calibrate upon installation
- Re-calibrate if does not pass verification or if instrument setup is modified

Verification

- Determine how often to verify data
- Run a reference material as a sample (in standard mode)
- Compare results vs literature values
- If results are within your tolerance – system checks out and does not need re-calibration
- If results are out of tolerance, then re-calibrate

41



41

Requirements Prior to Calibration

- DSC cell must be free of contaminants
- An inert purge gas, such as nitrogen, where the flow rate is controlled to 10-50 ml/min +/- 5 ml/min
- A balance to weigh specimens and containers to at least +/- 0.1 mg. The balance should have a capacity greater than 20 mg.
- High purity reference materials (>99.99%) for calibration

42



42

Instrument Setup Factors Affecting Calibration

Purge Gas Type

- Re-calibrate baseline/Tzero, temperature and cell constant
- Thermal conductivity of helium \neq Thermal conductivity of nitrogen/air/oxygen \neq Thermal conductivity of argon

Cooling Accessories

- Re-calibrate baseline/Tzero, temperature and cell constant
- The position of the cooling head around the cell will affect the calibration of the instrument. Uninstallation and reinstallation of a cooling accessory or changing the cooling accessory warrants a complete re-calibration

Pan Selection

- Re-calibrate temperature and cell constant
- It will not impact the baseline/Tzero calibration


43



43

Temperature and Cell Constant Calibration

- Prepare a 3-5 mg sample of indium and “pre-melt” prior to first use

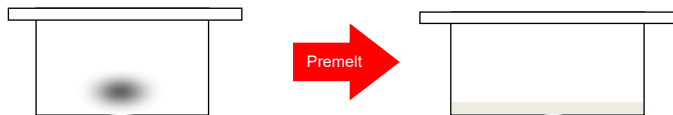

- Verify at least once a month
- Typical values for cell constant:
 - 0.9 to 1.2 (in N₂)
 - Helium will typically give higher values for the cell constant

44



44

Verifying Cell Constant & Temperature


- Run Indium as a sample (i.e. in standard mode not calibration mode)
- Analyze melt and record melt onset & heat of fusion
- Compare to known values
 - Melting of In 156.6°C
 - Heat of Fusion 28.71J/g

45



45

Verifying Cell Constant & Temperature

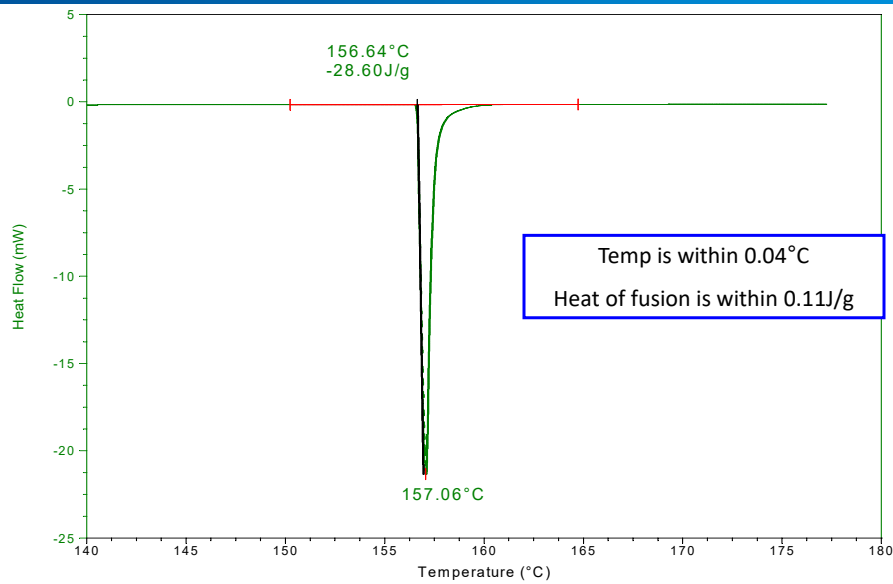

46



46

Reference Standards for Calibration

- Enthalpy (cell constant)
 - Benzoic acid (147.3 J/g) T_m = 123°C
 - Urea (241.8 J/g) T_m = 133°C
 - **Indium (28.71 J/g) T_m = 156.6°C**
 - Anthracene (161.9 J/g) T_m = 216°C
- Temperature
 - Cyclopentane* -150.77°C
 - Cyclopentane* -135.09°C
 - Cyclopentane* -93.43°C
 - Cyclohexane # -83°C
 - Gallium # 29.76°C
 - Phenyl Ether # 30°C
 - p-Nitrotoluene E 51.45°C
 - Naphthalene E 80.25°C
 - **Indium # 156.60°C**
 - Tin # 231.95°C
 - Lead* 327.46°C
 - Zinc # 419.53°C

* GEFTA recommended
Thermochim. Acta, 219 (1993) 333.

ITS 90 Fixed Point

E Zone refined organic compound
(sublimes)

47



47

Verifying Baseline

- Run Empty cell (no pans), -90°C to 400°C (w/ RCS) at 20°C/min
 - Experiment is run in the standard mode
 - Plot mW vs. temperature on a 1mW scale
 - Should look fairly flat on this scale
 - Should be around zero heat flow
 - Measure drift and compare to instrument specifications.
 - Verify performance periodically

48



48

Verifying Baseline


- Importance of a flat baseline:
 - Detecting very weak transitions
 - Accurate integration of enthalpy
 - Kinetics, partial area analysis, extent of reactions
 - Initial crystallinity

49



49

Verifying Baseline

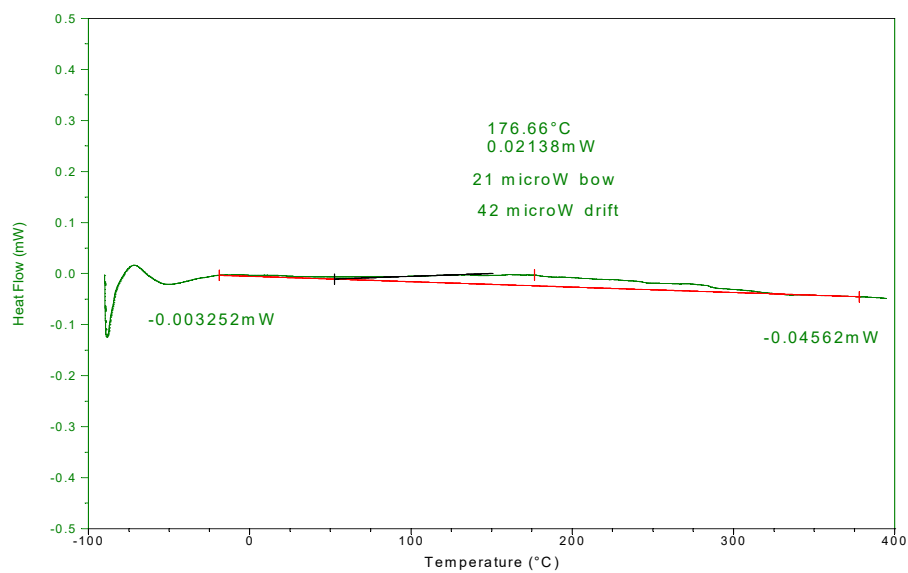


50



50

Empty Cell Baseline at 20 °C/min – DSC2500

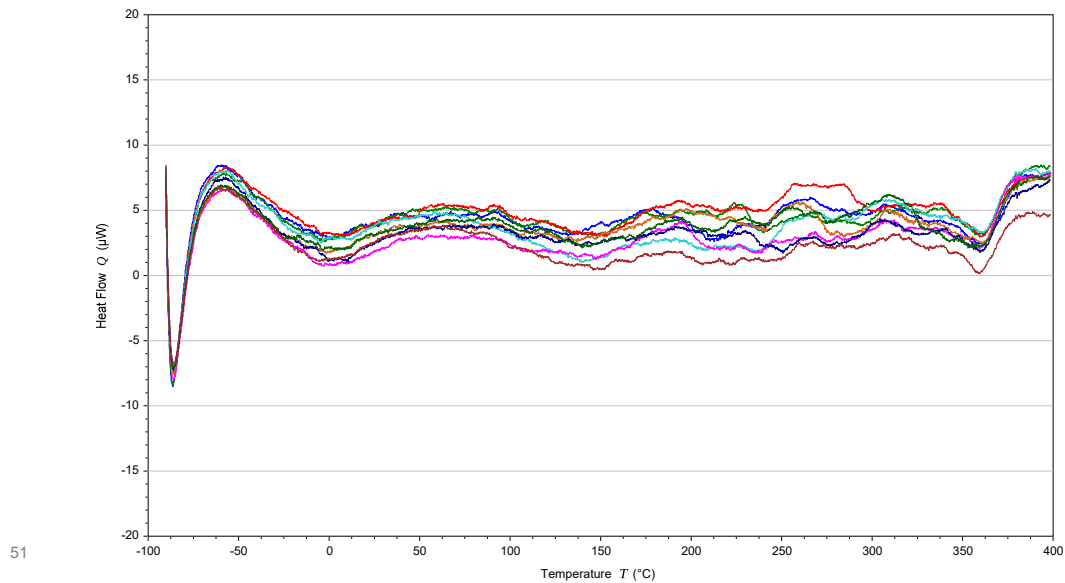

MDSC Calibration – Discovery Series


- Calibrate your DSC as normal
 - Tzero™
 - Cell Constant
 - Temperature
- Cp Calibration is Optional
- If measuring absolute quantitative Cp then...
- Need to calibrate Reversing Heat Capacity or Direct Heat Capacity (DSC2500 only)

52



52

Cp Calibration in Trios

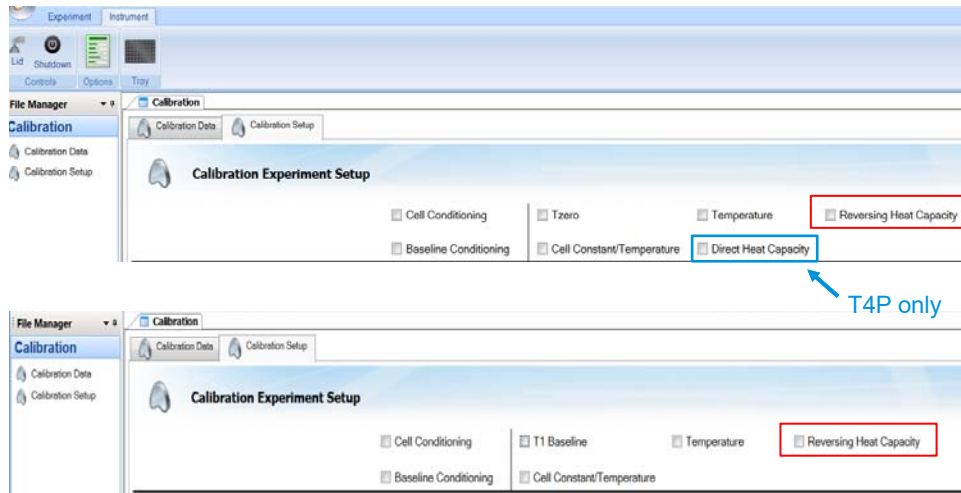


53



53

MDSC Cp calibration using TRIOS

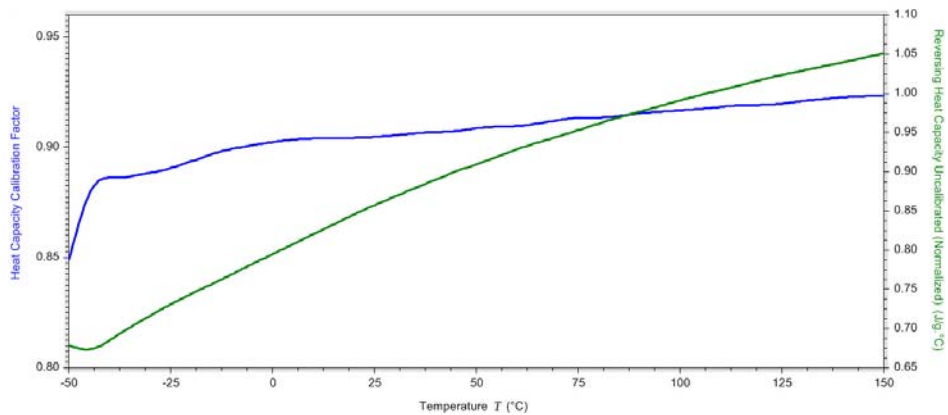


KCp is determined as a continuous function of temperature in the Discovery DSC.

54



54

MDSC C_p verification using TRIOS

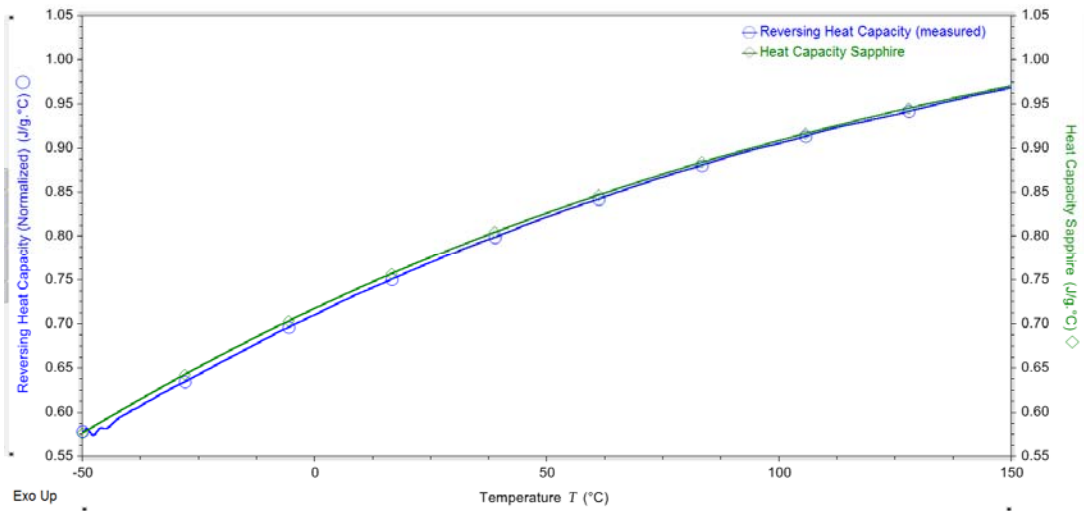

55



55

Experimental Design: Instrument Set Up



56

Instrument Hardware and Gas Selection Considerations



Temperature Range Dependent On The Cooling System

- Finned Air Cooling System (FACS): Ambient to 725°C
- Quench Cooling Accessory (QCA): -180°C to 400°C
- Liquid Nitrogen Cooling System (LN2P): -180°C to 550°C
- RCS120: -120°C to 400°C
- RCS90: -90°C to 550°C
- RCS40: -40°C to 400°C



Purge Gas Selection

- **Nitrogen**
 - inert, inexpensive and readily available
 - flow rate of 50ml/min
- **Helium**
 - a high thermal conductivity gas which improves response time and cooling capabilities
 - the recommended purge gas when using the LN2 accessory at temperatures below -100°C
 - flow rates of 10-25ml/min are typically used; cell constant affected by flow rate
- **Air/Oxygen**
 - used when studying oxidative stability of materials



Sample Press and Pan Selection

- Aluminum: max. temperature of 600°C
- Gold
- Copper
- Graphite, Alumina
- Platinum
- Stainless Steel


57



57

Cooling Accessories

- Finned Air Cooling System (FACS): Ambient to 725°C
- Quench Cooling Accessory (QCA): -180°C to 550°C *
- Liquid Nitrogen Cooling System (LNCS): -180°C to 550°C
- RCS120: -120 °C to 400 °C
- RCS90: -90°C to 550°C
- RCS40: -40°C to 400°C

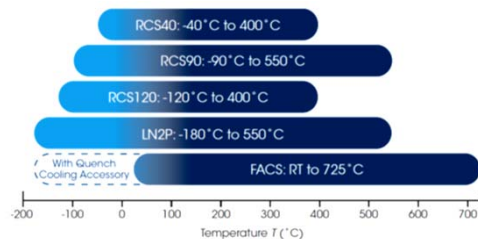

58

* Isothermal or slow heating ramp experiments above 400°C are not recommended.



58

Selecting the Cooler – Discovery DSC

TA Instruments TRIOS

Application
Discovery DSC
Information
General
Cooler
Auto Sampler
Temperature Cal
Heat Capacity

Cooler Settings

Cooler Selection This is used to select the cooler type

Cooler Selection: RCS 90 Cooler

Activate secondary purge when lid is opened (R)

Between Runs

Leave Cooler On

Autofill LN2P if below: 60.0 %

59

59

Selecting the Purge Gas – Discovery DSC

TA Instruments TRIOS

Experiment

Application
Discovery DSC
Information
General
Cooler
Auto Sampler
Temperature Cal
Heat Capacity

Global Settings

Global Options

Transition Direction: Exotherm Down Exotherm Up

Heat Flow Selection: Heat Flow T4P (mW)

Data Sampling Interval: 0.1 s/pt

Lid Type: Standard Temperature Lid

Enable Sequence without Using Autosampler

Gas Connections

Gas 1: Nitrogen

Gas 2: Nitrogen

Stop the experiment when Air

Instrument Database Backup (M): Nitrogen (04, 2015)

Backup: Daily Around 1:00 AM Backup Now

Backup Folder: C:\ProgramData\TA Instruments\TRIOS\ThermalDBBackups

Prohibit any scheduled backups on this computer

View Log Restore

Options Exit

60

This is used to specify the type of gas connected to Gas #1 and Gas #2 inlets

60

Setting the Purge Gas Flow Rate – Discovery DSC

This is used to select which gas is going to the DSC cell and the flow rate for that gas.

General

Temperature: 40.01 °C

Flange Temperature -82.61 °C

Cooler Selection RCS 90 Cooler

Standby 40.0 °C

RCS ON

Gas 1 Nitrogen: 50.0 mL/min

Base purge 263.4 mL/min

Gas 1: Nitrogen

Flow rate 50

61



61

Recommended Purge Gas Flow Rates

<u>Module</u>	<u>Purge Port</u>
All TA DSC's	50 ml/min (N ₂) or 25 ml/min (He)

- If purge gas is too slow - possible moisture accumulation & early aging of the cell
- If purge gas is too fast – excessive noise

62



62

Experimental design: Sample Preparation and Considerations



63

TGA for DSC Experimental Design

- Thermogravimetric Analysis (TGA) measures weight loss or gain as a function of temperature, time and atmosphere.
- General applications of TGA include:
 - thermal stability
 - residual solvent, out gassing, moisture sorption/desorption
 - filler/fiber content
 - weight loss on cure
- TGA measurements are extremely useful in selecting experimental conditions for DSC experiments and for interpreting results.




64



64

Typical TGA data: TGA of Drug A Monohydrate

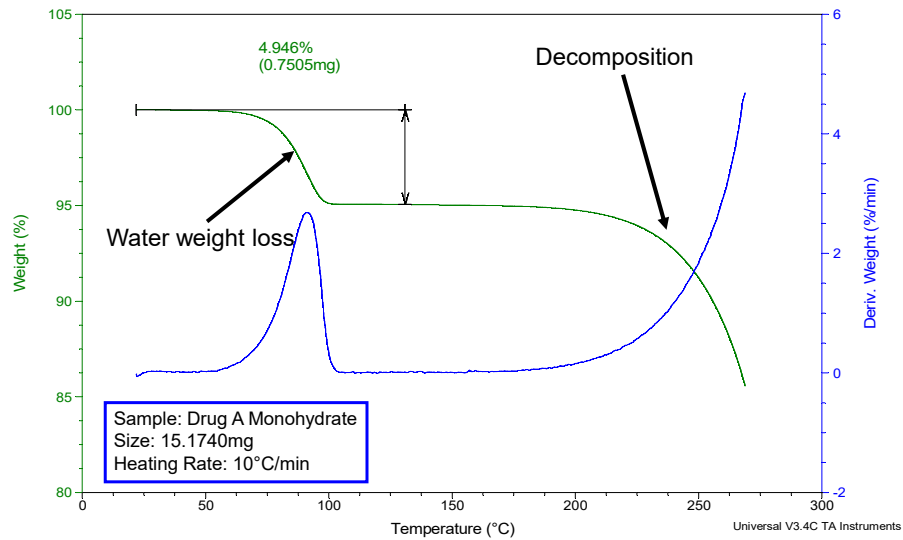

65



65

Selecting Optimum Experimental Conditions

- Use TGA data to help select DSC experimental conditions
 - Standard (non-hermetic) vs. Hermetic (sealed) pans
 - Use hermetic pan if sample loses approximately 0.5% weight or more.
 - Use hermetic pan with pin hole lids if sample loses volatiles such as water
 - Maximum Temperature
 - Excessive decomposition will contaminate the DSC cell between runs
 - When comparing samples, always use the same experimental conditions

66



66

DSC Pan Selection

- High thermal conductivity – aluminum, gold, copper, platinum
- Inert - alodined aluminum, ceramic, graphite
- Flatness of the pan for optimal thermal contact
- Standard, non-hermetic vs. hermetic sealing
- Capacity/sample volume
- Temperature range



67



67

Sample Pans

- Type of pan depends on:
 - Sample form
 - Volatilization
 - Temperature range
- Use lightest, flattest pan possible
- Always use reference pan of the same type as sample pan

68



68

Sample Shape

- Keep sample thin
- Cover as much the bottom of pan as possible



69



69

TA Instruments Tzero Pans (Aluminum)

Tzero Pan



Tzero Low-Mass Pan



- The Tzero pan has been engineered to have a perfectly flat bottom and not to deform during crimping. This ensures the optimal contact between pan and sensor, minimizing the contact resistance and improving resolution.
- The Tzero Pan can be configured for non-hermetic or hermetic use. P/N 901683.901 Tzero Pans (pkg. of 100)
- The Tzero Low-Mass Pan is designed for the highest sensitivity when sample mass is limited. P/N 901670.901 Tzero Low-Mass Pans (pkg. of 100). Can only be used with the non-hermetic Tzero lid.


70



70


TA Instruments Tzero Pans

Tzero Lid

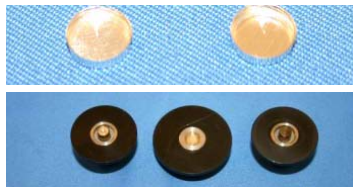

- Tzero Lid (P/N: 901671.901) - Lightweight aluminum lids for use in sample encapsulation with the Tzero Pans and the Tzero Low-Mass Pans. The seal is not hermetic.

Tzero Hermetic Lid

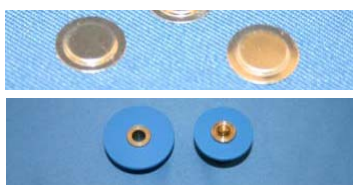

- Tzero Hermetic Lid (P/N: 901684.901) (pkg. of 100) and P/N: 901685.901 Tzero Hermetic Pinhole Lid (75 micron diameter pinhole) (pkg. of 50). Used only with the Tzero pan, not the low mass Tzero pan

71



71

Standard Series DSC Other Pans (Non Hermetic)


- Standard pans are available in:
 - Gold (p/n 900866.901 pan, p/n 900868.901 lid): up to 725°C
 - Graphite (p/n 900874.901 pan, p/n 900873.901 lid): up to 725°C (in N₂)
- Standard Pans with no lids available
 - Platinum (p/n 900578.901): up to 725°C
 - Copper (p/n 900867.901) : up to 725°C

72



72

Hermetic DSC Pans

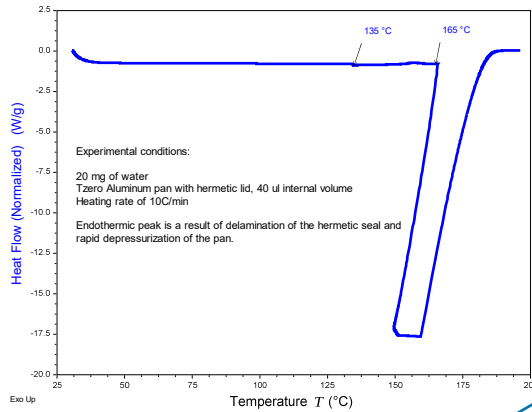

Aluminum
Hermetic
3 atm



Sample:
 1. Liquid
 2. Solid with
 volatile content

High Volume
Stainless
Steel
40 atm



Gold
Hermetic
6 atm



73





73

Hermetic Pans Summary

- Hermetic Pans are available in:
 - Aluminum: <600°C; <3 atm (300 kPa gauge)
 - Alodined Aluminum: <200°C; <3 atm (300 kPa gauge)
 - Gold: <725°C; <6 atm (600 kPa gauge)

- Specialized Sealed Pans

- High Volume: 100µL; <250°C; 37 atm (575 psi)
- P/N 900825.901
- High Pressure: 35µL; <300°C; 100 atm (1450 psi)
- P/N 900808.901

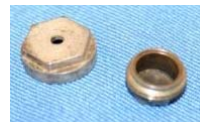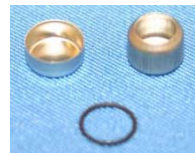


Note: 3 atm is approximately 44 psi

74



74

It Can Matter What Pan You Use

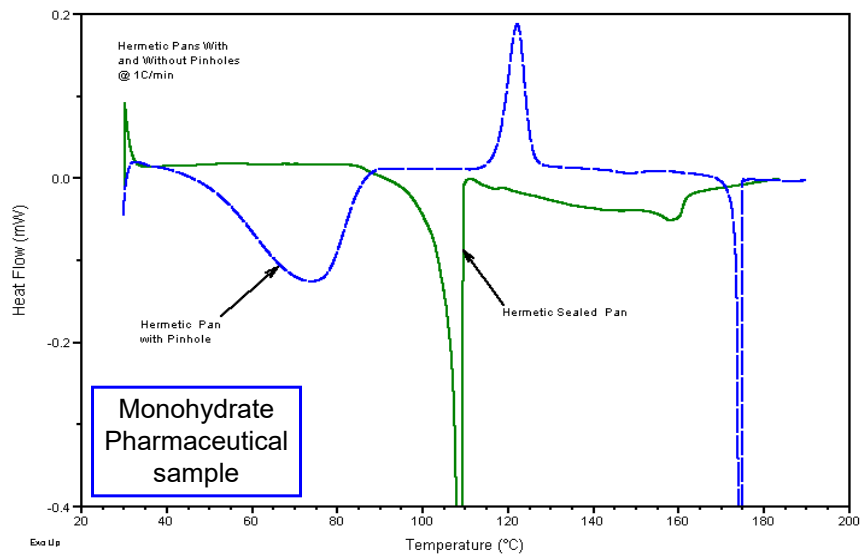

75

What if Sample Spills out of the Pan? Keeping the DSC Cell Clean

- One of the first steps to ensuring good data is to keep the DSC cell clean
- How do DSC cells get dirty?
 - Decomposing samples during DSC runs
 - Samples spilling out of the pan
 - Transfer from bottom of pan to sensor

76

TA

76

Cleaning the Cell

- Use solvent – slightly damp swab with an appropriate solvent
 - Heat cell to 200°C for 10 min to drive off any remaining solvent

- If the cell is still dirty
 - Clean w/ brush
 - Be careful with the Tzero™ thermocouple
 - Fibers in cell from cleaning brush need to be removed



77



77

Cleaning Cell: Bake Out Procedure

- Bake out
 - Should be used as a last resort if none of the previous steps are effective
 - Involves Air purge and/or an open lid
 - Heat at 20°C/min to appropriate temp (max of 550°C on Q series, max. 400°C in Discovery)
 - Do NOT hold Isothermal at the upper temperature
 - Cool back to room temp & brush cell again

- Irrespective of the cleaning method used, always verify the baseline at the end of the cleaning procedure, and recalibrate the DSC if required

- Check out the TA Tech tip video on cleaning the DSC cell:
 - <https://www.youtube.com/watch?v=cclJXrbUICA>

78



78

Method Development



79

DSC General Method Recommendations


- Run a Heat-Cool-Heat at 10-20 °C/min
- Use specific segments as needed, i.e. gas switch, abort, etc.
- Ensure that the starting temperature of the experiment is chosen to encompass the entire transition (see next slide)
- Modify heating rate based on what you're looking for

80



80

Have 2 Minutes of Baseline

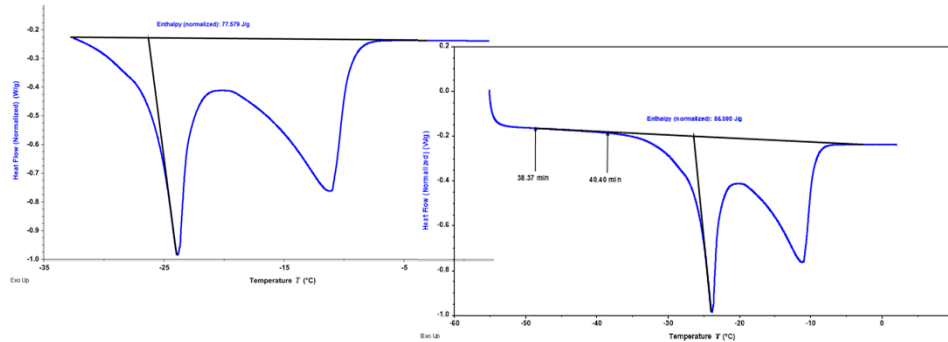

- Start Temperature
 - Generally, the baseline should have two (2) minutes to completely stabilize prior to the transition of interest. Therefore, at 10°C/min., start at least 20°C below the transition onset temperature
- End Temperature
 - Allow a two (2) minute baseline after the transition of interest in order to correctly select integration or analysis limits

81



81

Method Development

A Heat Cool Ramp Method


- 1) Ramp 10°C/min to -90°C or Equilibrate to -90°C
- 2) Ramp 10°C/min to 200°C
- 3) Ramp 10°C/min to -90°C
- 4) Ramp 10°C/min to 200°C

82



82

Heat Cool Heat Cycles of High Density Polyethylene (HDPE)

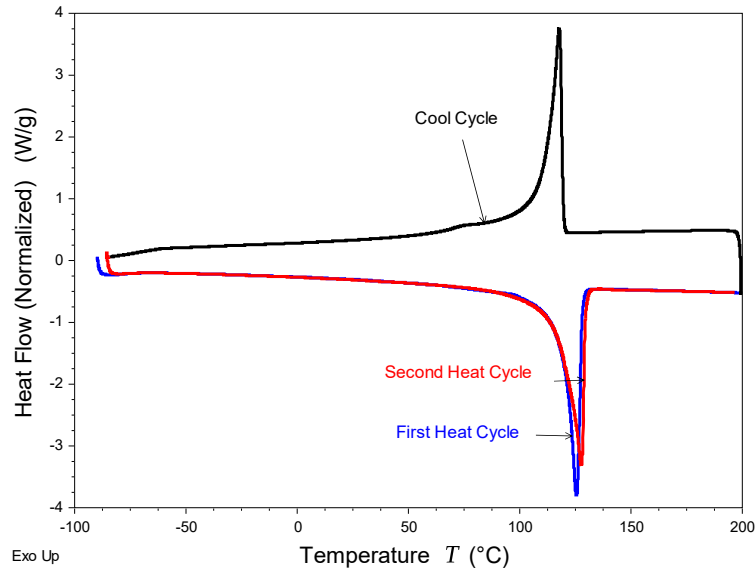

83



83

How to Program MDSC?

Period

- Typically 60 seconds for transitions
- Typically 120 seconds for Cp

Heating rate

- Typically 2-3°C/min

Amplitude

- Typically $\pm 1^\circ\text{C}$ for amorphous transitions and Cp
- Use Heat-Iso amplitude for crystalline transitions

84



84

Starting Guidelines for Determining the Tg in Amorphous Materials – Conventional MDSC®

- Conventional MDSC conditions
 - Period: 60 second
 - Minimum period depends on pan type and pan mass
 - Can use 40 sec for low-mass Tzero pans
 - Use 200 sec for Hi-volume pans
 - Amplitude: $\pm 1^\circ\text{C}$
 - Increase for added sensitivity
 - Heating Rate $3^\circ\text{C}/\text{min}$
 - Slow down heating rate for sharper transitions or large enthalpic recovery peaks

85



85

Starting Guidelines for crystallinity/melting studies: MDSC® Heat Only Conditions


- MDSC Heat-Only Conditions
 - Period: 60 second
 - Minimum period depends on pan type and pan mass
 - Can use 40 sec for low-mass Tzero pans
 - Heating Rate $2\text{-}3^\circ\text{C}/\text{min}$
 - Ensure 4-5 cycles through peak @ half-height of melt
 - Amplitude Heat-only (Heat-Only/Heat-Iso)

86



86

MDSC – Heat-only Temperature Modulation

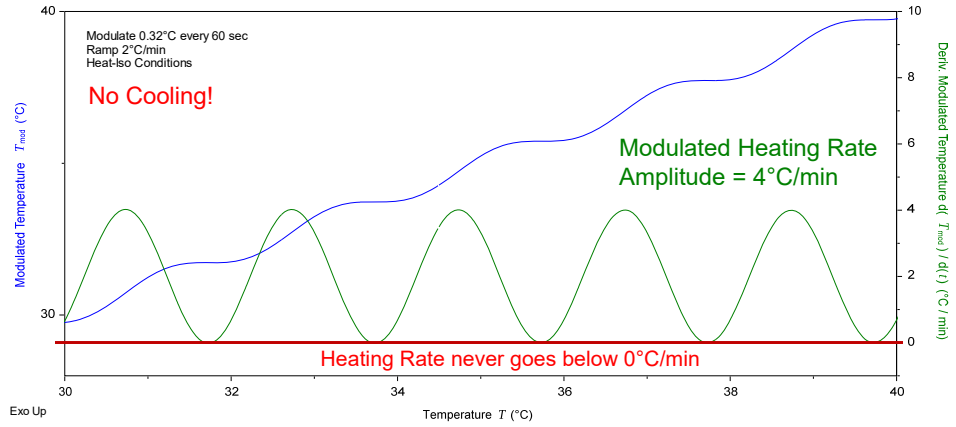


87



87

MDSC – Heat-Cool Temperature Modulation

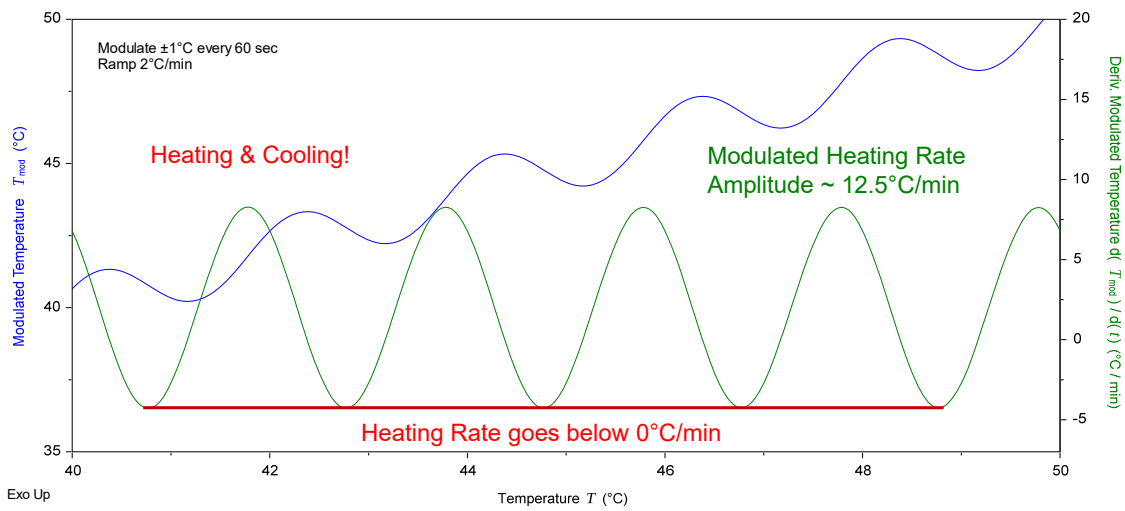

88



88

Parameters for Heat Capacity (Cp)

- **Conventional MDSC conditions**
 - Period 120 seconds
 - Can use 100 seconds with low-mass Tzero pans
 - Amplitude $\pm 1^\circ\text{C}$
 - Heating Rate $3^\circ\text{C}/\text{min}$
 - Can also do Quasi-isothermal for specific temperatures

89



89

MDSC Conditions for Discovery SDT 650

- Used to measure Heat Capacity
- Amplitude = $\pm 3^\circ\text{C}$
- Period = 200 seconds
- Heating Rate = $5^\circ\text{C}/\text{min}$
- Best above 300 or 400°C

90



90

MDSC Parameters: The Modulation Period

- The modulation period should be long enough to allow for a quantitative measurement of heat flow between the sensor and the sample.
 - Too short a period may result in poor data due to thermal resistances between the sensor and the pan, the pan and the sample, and the sample itself.
 - Too long a period requires a lower average heating rate resulting in potential signal to noise issues and longer experiment time.
- Suggested starting periods:
 - Thermal transitions: 60 seconds
 - Heat capacity measurement: 120 seconds
- The default conditions are also a good starting point for MDSC experiments.

91



91

Modulation Period: Summary Guidelines

- Choice of period will also be dependent to some extent on the pan you use:
- For transitions:
 - Use 40 sec for pans <30mg (Tzero Low Mass)
 - Use 60 sec for pans ~50mg (Tzero Std & Hermetic)
 - Use 200 sec for Hi-Volume Pans
- For quantitative Heat Capacity measurements:
 - Use 100 sec for pans <30mg (Tzero Low Mass)
 - Use 120 sec for pans ~50mg (Tzero Std & Hermetic)

92



92

MDSC[®] Parameters: The Modulation Amplitude

- Factors Considered Choosing Modulation Amplitude:
 - Signal to Noise: Larger amplitudes will result in greater sensitivity, but if too large may effectively smear a transition (loss of resolution) and may adversely affect the linearity of the heat flow response.
 - Cooling? – Amplitude is selected based on whether cooling (negative heating) is desired during the modulation.
 - Amplitudes utilizing cooling are effective for determining heat capacity or changes in heat capacity such as glass transitions.
 - Amplitudes resulting in only positive or minimum zero heating rates are effective in analysis of melt transitions.

93



93

MDSC Parameters: The Modulation Amplitude

- The software will allow choice of amplitudes between +/-0.001 and 10 °C.
- The practical range is +/- 0.1 and 2 °C, with lower amplitudes yielding poor sensitivity and high amplitudes decreasing resolution.

94



94

Modulation Amplitude: Summary Guidelines

- For amorphous transitions and measuring heat capacity:
 - Start with $\pm 1^\circ\text{C}$
 - Increase amplitude for added sensitivity
- For crystalline transitions:
 - Use Heat-Only

95



95

MDSC Parameters: Average Heating Rate

- The average heating rate is chosen to maximize signal to noise while balancing the need to obtain 4-6 modulations through the temperature range of interest.
- For melting transitions, the range of interest is determined at the half-height of the melting peak.
- For glass transitions, it is determined where the transition is changing fastest.
- Typically a heating rate from 1-3 $^\circ\text{C} / \text{min}$ is chosen.
- The heating rate is the best parameter to vary when 'fine tuning' the MDSC experiment.

96



96

Selecting MDSC[®] Underlying Heating Rate


- Check your data by plotting the modulated heat flow and the total heat flow signals. You should obtain 4-6 modulations through the transition of interest.
- Faster rates improve signal quality and sensitivity balanced by the need for sufficient modulations.
- Often MDSC optimization is an iterative process, so its not unusual to repeat an experiment.

97



97

MDSC Parameters: Sufficient Modulations Through Melting Transition

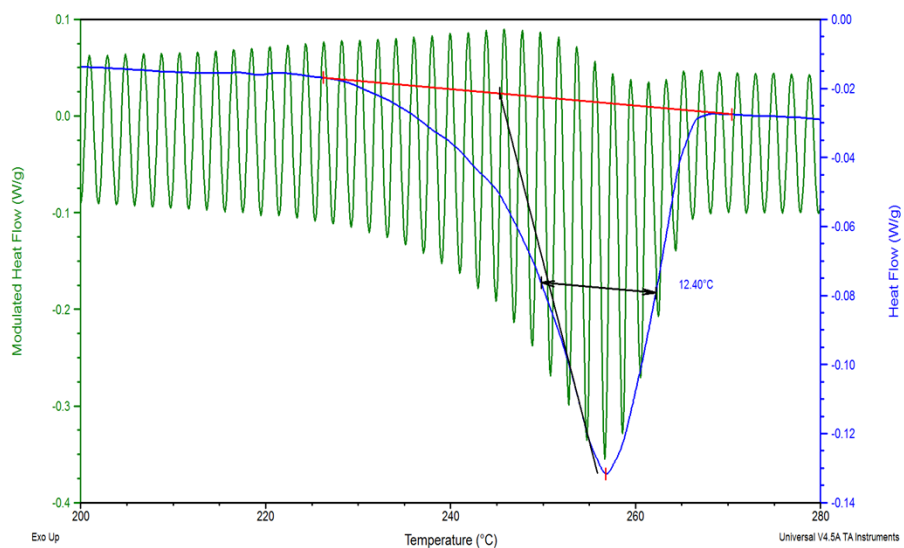


98



98

Selecting Signals for Plotting (TRIOS)

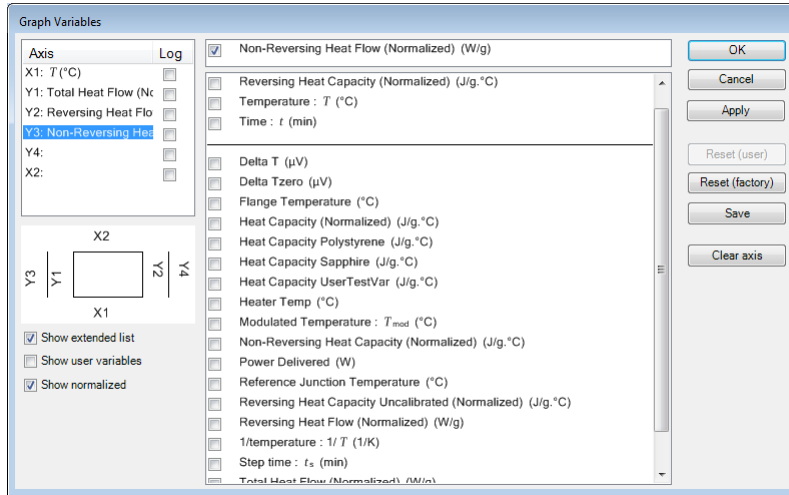


99



99

MDSC® of PET

Default Universal Analysis Presentation of Data

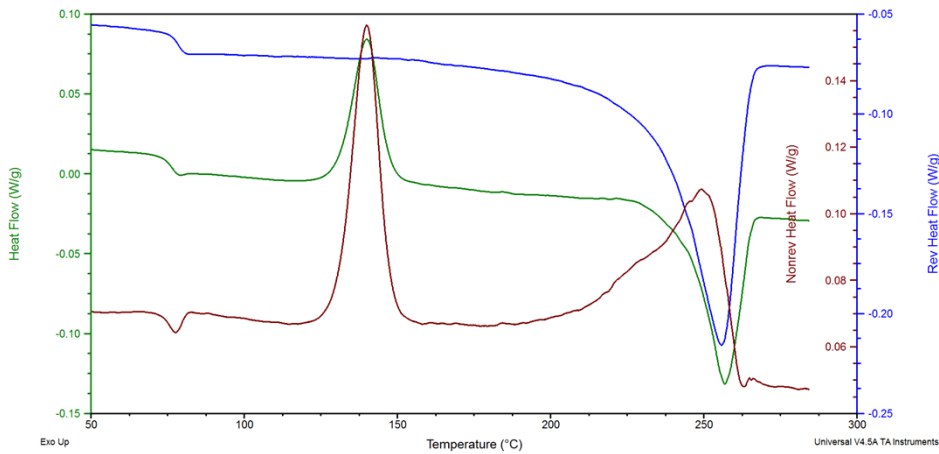


100



100

MDSC[®] of PET (Normalized Data)

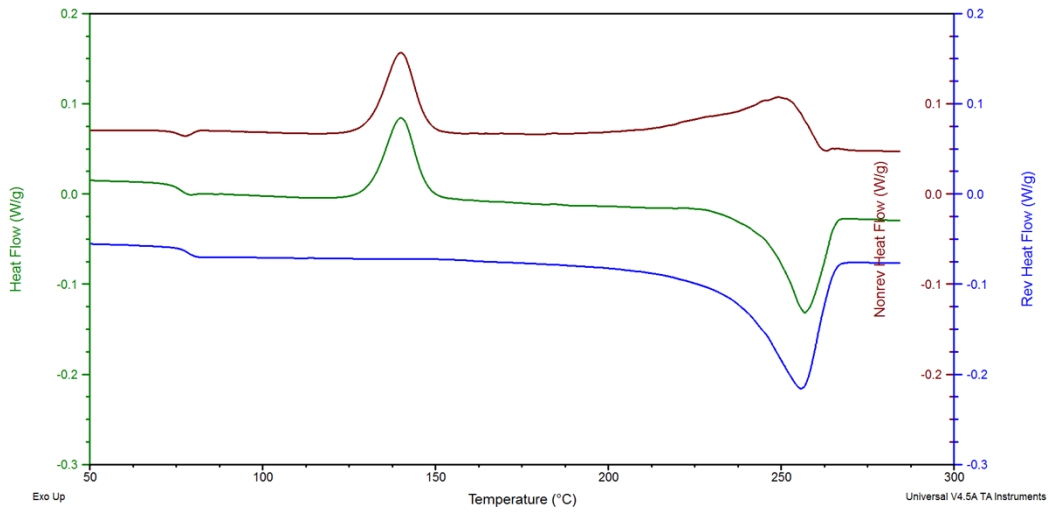


101



101

MDSC[®] of PET – Raw Signals

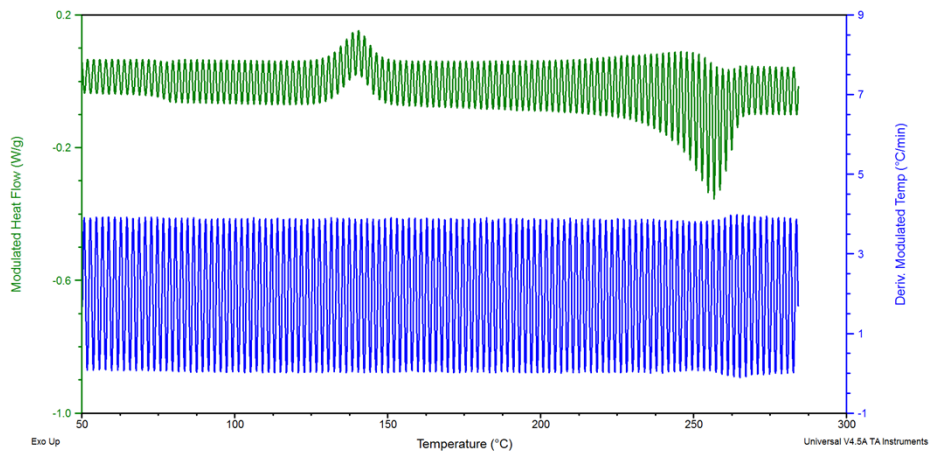


102



102

MDSC[®] of PET

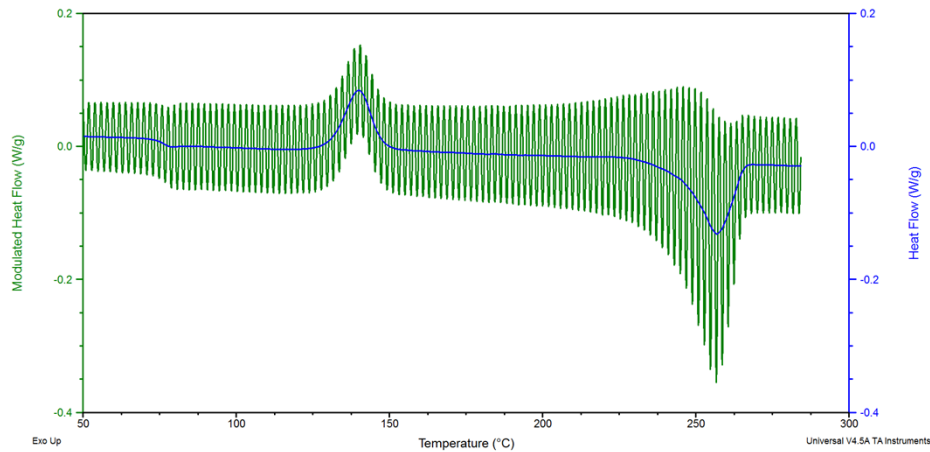

103



103

Summary for Method Development

- Typical starting sample mass – depends on the size of the transition
 - 5-10mg for most polymers
 - 3-5 for chemical melting
 - <1 for explosives
- Typical starting ramp rate = 10°C/min (Heat-Cool-Heat)
 - Faster for better sensitivity
 - Slower for better resolution
- Change above as needed
- MDSC[®] as needed

104



104

What if I need help?

- TA Tech Tips
 - <http://www.youtube.com/tatechtips>
- TA Instruments Applications Helpline available from the TA website
 - <http://www.tainstruments.com/support/applications/applications-hotline/>
- Check out our Website
 - <http://www.tainstruments.com/>

105



105

What if I need help?

- TA Tech Tips
 - <http://www.youtube.com/tatechtips>
- TA Instruments Applications Helpline available from the TA website
 - <http://www.tainstruments.com/support/applications/applications-hotline/>
- Check out our Website
 - <http://www.tainstruments.com/>

106



106

- For additional questions:
- Email thermalsupport@tainstruments.com
- Please put Online Training Questions in the subject line
- You can download this presentation from:
 - <https://www.tainstruments.com/online-training-course-downloads/>



107

TA Instruments –Waters LLC
www.tainstruments.com



108